(DE) (AT)	Heizkissen	(ES)	Almohadilla eléctrica
(CH)	<i>Gebrauchsanweisung</i>2		<i>Instrucciones para el uso</i>43
(FR) (BE)	Coussin chauffant	(DK)	Varmepude
(CH)	<i>Mode d'emploi</i>9		<i>Betjeningsvejledning</i>50
(NL) (BE)	Warmtekussen	(HU)	Melegítő párna
	<i>Gebruiksaanwijzing</i>16		<i>Használati útmutató</i>56
(PL)	Poduszka elektryczna	(IT) (MT)	Termoforo
	<i>Instrukcja obsługi</i>23	(CH)	<i>Istruzioni per l'uso</i>64
(CZ)	Výhřevná poduška	(GB) (MT)	Heat pad
	<i>Návod k použití</i>30		<i>Instruction for Use</i>71
(SK)	Elektrická poduška		
	<i>Návod na obsluhu</i>36		

Service-Hotline:

(DE) Tel.: 0800 724 2355
(AT) Tel.: 0800 212 288
(CH) Tel.: 0800 200 510
(BE) Tél.: 0800 706 11
(NL) Tel.: 0800 543 0543

(FR) Tél.: 0805 340 006
(PL) Tel.: 800 707 009
(CZ) Tel.: 0800 555 013
(SK) Tel.: 0800 606 018
(ES) Teléfono: 800 009 222

(DK) Tlf.: 8070 5789
(HU) Tel.: 0680 106 886
(IT) Tel.: 800 141 010
(GB) Tel.: 0800 931 0319
(MT) Tel.: 800 626 05

Inhalt

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Contents

1. Items included in the package	72	4.4 Setting the temperature.....	74
1.1 Description	72	4.5 Automatic switch-off	75
2. Important instructions		4.6 Switching off	75
Retain for future use	72	5. Cleaning and maintenance	75
3. Intended use	74	6. Storage	75
4. Operation	74	7. Disposal	76
4.1 Safety	74	8. What if there are problems?	76
4.2 Initial use	74	9. Technical data	76
4.3 Switching on	74	10. Guarantee/service	76

Explanation of symbols

The following symbols are used on the device, in these instructions for use, on the packaging and on the type plate for the device:

	WARNING: Warning of risks of injury or health hazards		
	CAUTION: Safety information about possible damage to appliance/accessories		
	NOTE: Important information		
	Read the instructions		Do not insert pins
	Do not use folded or rucked		Not to be used by very young children (0 - 3 years)
	Manufacturer		
	The device has double protective insulation and therefore complies with protection class 2.		
	Please dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment).		
	Marking to identify the packaging material. A = Material code, B = Material number: 1-7 = Plastics, 20-22 = Paper and cardboard		
	This product satisfies the requirements of the applicable European and national directives.		
	Manufacturer code for production data		
	Separate the product and packaging elements and dispose of them in accordance with local regulations.		
	This device meets the high requirements of OEKO-TEX® STANDARD 100 for human ecology as verified by the Hohenstein Research Institute.		


	Wash at a maximum temperature of 30 °C, Very gentle wash		Do not iron
			Do not bleach
	Do not dry-clean		Do not dry in a tumble dryer

1. Items included in the package

- 1 Heat pad
- 1 Controller
- 1 Instructions for use

1.1 Description

- 1 Mains plug
- 2 Mains cable
- 3 Controller
- 4 Indicator lamp
- 5 Sliding switch for ON/OFF and temperature settings
- 6 Plug-in coupling

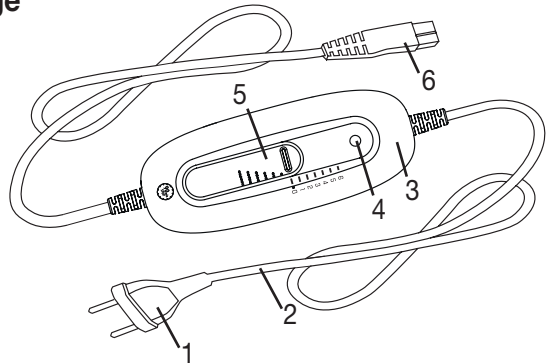2. Important instructions

Retain for future use

WARNING

- Non-observance of the following notes may result in personal injury or material damage (electric shock, skin burns, fire). The following safety and hazard information is not only intended to protect your health and the health of others, it should also protect the product. For this reason, pay attention to these safety notes and include these instructions when handing over the product to others.
- This heat pad must not be used by people who are not sensitive to heat or by vulnerable people who may not be able to react to overheating (e.g. diabetics, persons with skin changes related to illness or persons with scarred skin in the area where the pad is to be applied, after taking any pain-relieving medication or alcohol).
- This heat pad may not be used by very young children (0-3 years old) as they are unable to respond to overheating.
- The heat pad can be used by children older than 3 and younger than 8 years of age provided they are supervised. For this, the control must always be set to the minimum temperature.
- This heat pad may be used by children over the age of 8 and by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the heat pad safely, and are fully aware of the consequent risks of use.

- If the heat pad is being used over several hours, we recommend that you set the lowest temperature setting on the control to avoid overheating the heated body part, which may lead to burns to the skin.
- Children must not play with the heat pad.
- Cleaning and user maintenance must not be performed by children unless supervised.
- This heat pad is not designed for use in hospitals.
- This heat pad is only intended for domestic/private use, not for commercial use.
- Do not insert needles.
- Do not use when folded or bunched up.
- Do not use if wet.
- This heat pad may only be used in conjunction with the controller specified on the label.
- This heat pad must only be connected to the mains voltage that is specified on the label.
- The electrical and magnetic fields emitted by this heat pad may interfere with the function of a pacemaker. However, they are still well within the limits: electrical field strength: max. 5000 V/m, magnetic field strength: max. 80 A/m, magnetic flux density: max. 0.1 millitesla.

Please consult your doctor and the manufacturer of your pacemaker before using this heat pad.

- Do not pull, twist or make sharp bends in the cables.
- If they are not positioned properly, there may be a risk of becoming entangled in, being strangled by, tripping over, or stepping on the cable and control of the heat pad. The user must ensure that excess lengths of cable, and cables in general, are safely routed.
- Please check this heat pad frequently for signs of wear and tear or damage. If any such signs are evident, if the heat pad has been used incorrectly or if it no longer heats up, it must be checked by the manufacturer before being switched on again.
- Under no circumstances should you open or repair the heat pad (including the accessories) yourself because faultless functionality can no longer be guaranteed thereafter. Failure to observe this will invalidate the warranty.
- If the mains connection cable of this heat pad is damaged, it must be disposed of. If it cannot be removed, the heat pad must be disposed of.
- Do not expose control and cables to direct sunlight. Otherwise the heat pad may get damaged.

- While the heat pad is switched on, the following conditions must be observed:
 - Do not place any sharp objects on it
 - Do not place any heat sources, such as hot water bottles, heat pads or similar, on it.
- The electronic components in the controller warm up when the heat pad is in use. For this reason, the controller must never be covered or placed on the heat pad when it is in use.
- Allow the heat pad to cool down before storing it away. Otherwise the heat pad may get damaged.
- To avoid any creasing in the heat pad, do not place any objects on top of it while it is being stored. Otherwise the heat pad may get damaged.
- It is essential to observe the information relating to the following chapters: Operation, Cleaning and maintenance, Storage.
- If you should have any questions about using our devices, please contact Customer Services.

3. Intended use

CAUTION

This heat pad is only designed to warm up the human body.

4. Operation

4.1 Safety

CAUTION

- The heat pad is fitted with a SAFETY SYSTEM. This sensor technology provides protection against overheating across the entire surface of the heat pad with an automatic switch-off in the event of a fault. If the SAFETY SYSTEM has switched off the heat pad, the indicator lamp (4) is no longer illuminated, even when the heat pad is switched on.
- Please note that for safety reasons, the heat pad can no longer be operated after a fault has occurred and must be sent to the specified service address.
- Do not connect the defective heat pad with another controller of the same type. This would trigger a permanent switch-off via the controller safety system.


4.2 Initial use

CAUTION

Ensure that the heat pad will not bunch up or become folded during use.

- First, connect the controller to the heat pad by inserting the plug-in coupling.
- Then plug the mains plug into the socket.

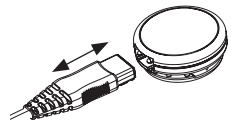4.3 Switching on

To switch the heat pad on, set the sliding switch for ON/OFF and temperature settings to setting 1, 2, 3, 4, 5 or 6. When it is switched on, the indicator lamp (4) lights up.

4.4 Setting the temperature

- Level 0: OFF
- Level 1: minimum heat
- Level 2-5: individual heat
- Level 6: maximum heat

**NOTE**

The fastest way to warm up the heat pad is to initially set the highest temperature setting.

**NOTE**

These heat pads have a fast heating function, which allows the pad to warm up faster in the first 10 minutes.

**WARNING**

If the heat pad is being used over several hours, we recommend that you set the lowest temperature setting on the controller to avoid overheating the heated body part, which may lead to burns to the skin.

4.5 Automatic switch-off

This heat pad is equipped with an automatic switch-off function. This turns off the heat supply approx. 90 minutes after the initial use of the heat pad. The indicator lamp (4) then begins to flash. So that the heat pad can be switched back on, the sliding switch for ON/OFF and temperature settings (5) must first be set to setting "0" (OFF). After about 5 seconds it is possible to switch it on again.

**4.6 Switching off**

To switch the heat pad off, set the sliding switch for ON/OFF and temperature settings to setting "0" (OFF). When it is switched off, the indicator lamp (4) no longer lights up.

**NOTE:**

If the heat pad is not in use, switch the sliding switch for ON/OFF and temperature settings to setting "0" (OFF) and unplug the mains plug from the socket. Then disconnect the controller from the heat pad by unplugging the plug-in coupling.

5. Cleaning and maintenance**WARNING**

Before cleaning, always remove the mains plug from the socket first. Then disconnect the controller from the heat pad by unplugging the plug-in coupling. Otherwise there is a risk of electric shock.

**CAUTION**

The controller must never come into contact with water or other liquids, as this may cause damage.

- To clean the controller, use a dry, lint-free cloth. Do not use any chemical or abrasive cleaning agents.
- Small marks on the heat pad can be removed with a damp cloth and if necessary, with a little liquid detergent for delicate laundry.

**CAUTION**

Please note that the heat pad may not be chemically cleaned, wrung out, tumble dried, put through a mangle or ironed. Otherwise the heat pad may be damaged.

- This heat pad is machine-washable.
- Set the washing machine to a particularly gentle wash cycle at 30 °C (wool cycle). Use a delicate laundry detergent and measure it out according to the manufacturer's instructions.
- Immediately after washing, reshape the heat pad to its original dimensions while it is still damp and spread it out flat on a clothes horse to dry.

**CAUTION**

- Do not use pegs or similar items to attach the heat pad to the clothes horse. Otherwise the heat pad may be damaged.
- Do not reconnect the switch to the heat pad until the plug-in connection and the heat pad are completely dry. Otherwise the heat pad may be damaged.

**WARNING**

Never switch the heat pad on to dry it! Otherwise there is a risk of electric shock.

6. Storage

If you do not plan to use the heat pad for a long period of time, we recommend that you store it in the original packaging. For this purpose, disconnect the controller from the heat pad by unplugging the plug-in coupling.

7. Disposal

For environmental reasons, do not dispose of the device in the household waste at the end of its useful life. Dispose of the unit at a suitable local collection or recycling point. Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal.

8. What if there are problems?

Problem	Cause	Solution
The indicator lamp is not illuminated while - the controller is properly connected to the heat pad - the mains plug is connected to a working socket - the sliding switch is set to 1, 2, 3, 4, 5 or 6	The safety system has switched the heat pad off permanently.	Send the heat pad and the controller for servicing.

9. Technical data

See the rating label on the heat pad.

10. Guarantee/service

Hans Dinslage GmbH, Riedlinger Straße 28, 88524 Uttenweiler (hereinafter referred to as “HaDi”), provides a guarantee for this product, subject to the requirements below and to the extent described as follows.

Claiming the guarantee shall not affect the Purchaser’s statutory guarantee rights arising from the purchase agreement with the Seller in the event of defects. The Purchaser may exercise these statutory guarantee rights free of charge. The guarantee shall apply without prejudice to any mandatory statutory provisions on liability.

HaDi guarantees the perfect functionality and completeness of this product.

The worldwide guarantee period is 3 years, commencing from the purchase of the new, unused product from the seller.

The guarantee only applies to products purchased by the buyer as a consumer and used exclusively for personal purposes in the context of domestic use. German law shall apply.

During the guarantee period, should this product prove to be incomplete or defective in functionality in accordance with the following provisions, HaDi shall carry out a repair or a replacement delivery free of charge, in accordance with these guarantee conditions.

Should the buyer wish to file a guarantee claim, they shall first contact the HaDi Customer Service team:

Service hotline (free):

☎ Tel.: 0800 931 0319

☎ Tel.: 800 626 05

E-mail:

service-uk@sanitas-online.de

service-mt@sanitas-online.de

The buyer will then receive further information about the processing of the guarantee claim, e.g. where they can send the product and what documentation is required.

If we request the buyer to send the defective product, the product must be sent to the following address:

NU Service GmbH, Lessingstraße 10 b, 89231 Neu-Ulm, Germany

A guarantee claim shall only be considered if the buyer can provide

HaDi, or an authorised HaDi partner, with

- a copy of the invoice/purchase receipt, and
- the original product.

The following are explicitly excluded from this guarantee:

- deterioration due to normal use or consumption of the product;
- accessories supplied with this product which are worn out or used up through proper use (e.g. batteries, rechargeable batteries, cuffs, seals, electrodes, light sources, attachments and nebuliser accessories);
- products that are used, cleaned, stored or maintained improperly and/or contrary to the provisions of the instructions for use, as well as products that have been opened, repaired or modified by the buyer or by a service centre not authorised by HaDi;
- damage that arises during transport between manufacturer and customer, or between service centre and customer;
- products purchased as seconds or as used goods;
- consequential damage arising from a fault in this product (however, in this case, claims may exist arising from product liability or other compulsory statutory liability provisions).

Repairs or an exchange in full do not extend the guarantee period under any circumstances.